



Laerdal[®]

helping save lives



When every second counts...

Thomas[®] Select[™] Tube Holder

In the initial critical moments of securing a patient's airway you need a device that you can rely on. The Thomas[®] Select[™] Tube Holder is designed to help you secure any airway device and ensure continued management of the patient's airway whether you use an endotracheal tube or supraglottic device. With the Thomas Select you now have the one go-to device to ensure this critical step of securing the airway is complete.

Quick and easy to use, the Thomas Select Tube Holder features a large opening combined with a push-screw mechanism specifically designed to secure a wide range of advanced airways, enabling you to spend less time securing the airway device and more time on what counts – saving lives.




www.laerdal.com


Securing SGA's (Supraglottic Airways) is recommended by the Journal of Emergency Medical Services as a key component for safe pre-hospital intubation

Thomas® Select™ Tube Holder


New features

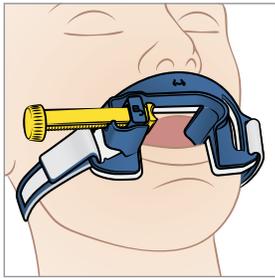


1. Larger opening
The Thomas Select Tube Holder is similar in shape to the Thomas Tube Holder, but has a larger opening to accommodate the wide range of airway tubes on the market. Includes position indicator symbol for correct placement.

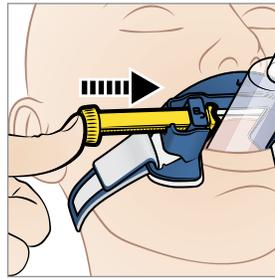

2. Push-screw mechanism
Quick and easy to use, a push-screw mechanism on the device for main adjustment. It is no longer necessary to screw all the way to tighten the tube, saving you time when you most need it.

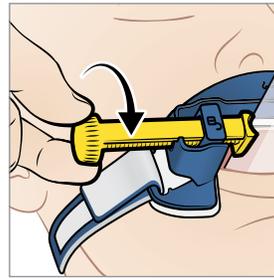

3. Normal screw mechanism
For final securement of the tube, the screw can be tightened as usual.

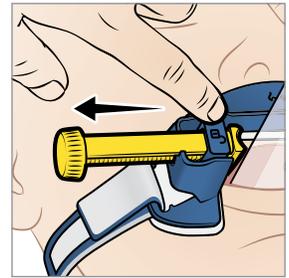

4. Quick Release
For a quick release of the tube clamp, simply push the unlock button.

Additional features and benefits

Large mouthpiece surface area	<ul style="list-style-type: none"> - Accommodates patients with removed dentures or facial trauma - Large opening allows for view of patient's mouth and lip color - Easy access for suctioning and other oral care needs
Wide Head Strap	<ul style="list-style-type: none"> - Provides secure and comfortable fit; even when wet
Pliable Bite Block	<ul style="list-style-type: none"> - Unique construction to help prevent the patient from biting or pinching the tube and restricting airflow. Pliable material to protect the teeth and gums.

Supplemental information

Size range for supraglottic airway devices: Size 3-6
Size range for tracheal tubes: ≥ 6.5 mm inner diameter
Accommodates tube sizes from > 6.5 mm (ID) to 42 mm (OD).

- Designed for single patient use only
- Not made with natural rubber latex
- Not to be used on children less than 8 years old

Ordering information:
600-42500 Thomas Select Tube Holder Adult, 25 pkg



For more information, visit www.laerdal.com



Laerdal[®]
helping save lives